NSK

Cleaning Adaptor

CLA-10

OPERATION MANUAL

For Varios2 Handpiece

CE Rx Only MADE IN JAPAN

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EN In this operation manual, VA2-LUX-HP Ti, VA2-HP and VA2-LUX-HP, are called "ultrasonic handpiece".
See the "Symbol" section on the back page for reference.

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EN User: Personnel with knowledge and skills required for operating this product.
Intended purpose: The ultrasonic handpiece cleaning adaptor CLA-10 is intended to clean ultrasonic handpieces in combination with the following products:
CLA-10 can be mounted to the washer-disinfector gasket (*1) and ultrasonic handpiece of Ultrasonic System manufactured by Nakanishi.

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EN **Related manuals**
For proper use with the combined items, refer to the following operation manuals:

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EN **Applicable product**
Products that can be used in combination with the device.

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⚠ PRECAUTIONS

EN

- Read this Operation Manual before use to fully understand the product functions and file for future reference.
- Do not use this cleaning adaptor for anything other than its intended purpose.
- Use the cleaning adaptor according to the intended purpose and operation procedure described in the operation manual by qualified personnel.
- Users are responsible for the operation control, maintenance and continual inspections of the cleaning adaptor.
- Do not drop, hit, or subject to excessive shock. The cleaning adaptor may be deformed.
- Do not disassemble or alter the cleaning adaptor.
- Be sure to check the cleaning adaptor for cleanliness and damage before use.
- Only NSK ultrasonic handpiece should be attached to the cleaning adaptor. Do not attach any other device.
- The cleaning adaptor may be damaged or discolored due to the effect of heat during mechanical cleaning. Replace with a new one immediately.
- This cleaning adaptor can be reused according to this operation manual.
- The O-ring is a consumable item and can be used for 3 months or 250 times. This cleaning adaptor (excluding the O-ring) is a consumable item and can be used 1,000 times. However, damage, discoloration, and deformation of the product may occur depending on the use environment.When the part with the slit that holds the handpiece becomes deformed, the holding force may become weak. If such a condition is found, replace the product with a new one. Perform periodic maintenance check. Repairs cannot be made by the user.
- Do not use the following liquids to wipe, immerse, or clean the product: strong/super acid water, strong acid/alkaline chemicals, chlorine-containing solutions, and solvents such as benzine or thinner.
- To ensure proper function and safe use of the medical device, perform periodic maintenance checks in accordance with the operation manual.
- If the washer-disinfector gasket (*1) provided by the user is worn out or deformed, the cleaning adaptor may come off. Periodically check the washer-disinfector gasket and the cleaning adaptor. If the washer-disinfector gasket is deformed or the cleaning adaptor comes off during cleaning, replace the washer-disinfector gasket with a new one.
- For reprocessing method of the handpiece, refer to OM-DE0981MA / OM-DE1011EN Varios2 Ultrasonic Handpiece OPERATION MANUAL.
- No special training is required for this product.
- Use appropriate disinfectant alcohol according to the regulations and guidelines of each country.
- Use only the cleaning and disinfection agent products specified in this operation manual.
- Wear eye protection, protective mask, and gloves for safety and infection prevention before starting the maintenance.
- Please report any serious incident that may occur in connection with this product and its accessories to the manufacturer as well as to your national authority.
- U.S. Federal law restricts this device to sale by or on the order of a licensed physician.

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EN

- Insert the cleaning adaptor into the washer-disinfector gasket until it stops (Fig. 1).
- When removing the cleaning adaptor from the washer-disinfector gasket, pull it upward while twisting it (Fig. 2).

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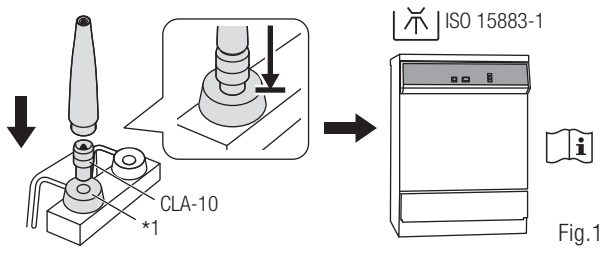
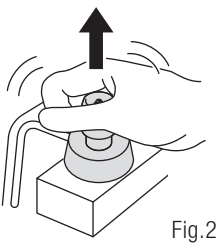
Fig.1

Fig.2

EN For the conditions for using the washer-disinfector, refer to Varios2 Ultrasonic Handpiece OPERATION MANUAL.

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PRECAUTIONS

EN
- After the thermo disinfection, take the product out of the washer-disinfector immediately (within an hour). Failure to do so causes corrosion.
- After using the washer-disinfector, dry the product completely, because remaining moisture may cause internal corrosion, etc.
- Before removing the cleaning adaptor from the washer-disinfector, make sure that the cleaning adaptor has cooled down sufficiently.

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EN During long-term use of the cleaning adaptor in a washer disinfector, the cleaning solution will stick to the surface of the cleaning adaptor.
Periodically remove the adaptor from the washer disinfector and clean the surface of the adaptor using a lint-free cloth dampened with disinfectant.
This product can be cleaned 1,000 times, dependeing on the use environment. However, if the product is damaged or deformed, replace it with new one.
Make sure that there is no deformation, damage or discoloration every three months. If any abnormalities are found, replace with a new one immediately.

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NSK Europe GmbH EC REP
NSK America Corp. www.nskdental.com