
NSK

OM-T0517E 009

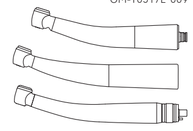OPERATION MANUAL

Ti-Max Z
Ti-Max X450
S-Max M
S-Max pico

Pana-Max2
Pana-Max PLUS
Pana-Max PLUS 2
FX Plus

MADE IN JAPAN

Download operation manuals

www.nsk-dental.com/support/userguide.html

NAKANISHI INC.
700 Shimohinata, Kanuma, Tochigi 322-8666, Japan
www.nsk-dental.com

NSK Europe GmbH
Elly-Beinhorn-Str. 8, 65760 Eschborn, Germany

EN

- User: Qualified professionals (dentist) • Intended purpose: The dental handpiece is intended for the following application(s): Caries removal, cavity and crown preparation, removal of dental restorations (fillings and prostheses), and finishing of teeth and dental restorations (preparation and adjustment), Only for Ti-Max X450 series : The Ti-Max X450 series is pneumatic dental drilling system handpiece with intended use of being a surgical tool for impacted third molar removal and periodontal procedures for which a conventional handpiece would be used. • Patient population : Any patient (ages 2 and over), The user's discretion is required prior to use.

- Depressing the Push Button while the handpiece is in rotation may lead to overheating, causing burn injuries or product failure, Avoid the push button to contact with any oral tissue.

- When operating the product always consider the safety of the patient. • Users are responsible for the operational control, maintenance and continual inspection of this product. • Do not disassemble or alter the handpiece and cartridge, etc. except as recommended by NSK in this Operation Manual. • Do not drop the handpiece from a height of more than 1 m. • Should the handpiece function abnormally during use, stop using the handpiece immediately and contact your Authorized NSK Dealer. • The product is delivered in a non-sterile conditions.Be sure to sterilize the product and the bur prior to use. • Follow the specifications described in this operation manual as to the bur to be used, etc. to prevent accidental release or breakage of the bur during rotation. • Perform regular function and maintenance checks. • To avoid clinical downtime it is recommended that a spare be kept on hand in case of a breakdown during surgery. • Avoid continual eye contact with the product. • Contact your Authorized NSK Dealer if repairs are necessary. • An accidental symptom, such as subcutaneous emphysema may be developed in the use of this product for extraction of wisdom teeth.Do not direct air exhausted from the handpiece to the bleeding area (wound) or into periodontal pockets. • Report all serious events related to the product immediately to the manufacturer and the responsible national authority. • Do not apply excess pressure to the bur as it may break or bend and become difficult to remove. • Always follow the instructions provided by the bur manufacturer regarding the bur usage and rotation speed, etc.

- This product is equipped with the Quick Stop Function, which means that when the air supply is stopped, the rotation stops in a shorter time than conventional handpieces, Also, when you turn the attached bur by hand, there will be resistance, but this is not a malfunction. • This product has a structure that a small amount of air escapes from the joint with the coupling, If liquid such as lubricating oil remains inside the handpiece, the liquid may leak from the joint, but this is not a malfunction. (S900, Z900L / Z890L / Z890KL / Z900KL / Z900SL) - Substances which are carcinogenic, mutagenic or toxic to reproduction (CMR), of category 1A or 1B, in accordance with Part 3 of Annex VI to Regulation (EC) No 1272/2008 of the European Parliament and of the Council, or substances having endocrine-disrupting properties for which there is scientific evidence of probable serious effects to human health and which are identified either in accordance with the procedure set out in Article 59 of Regulation (EC) No 1907/2006 of the European Parliament and of the Council or, once a delegated act has been adopted by the Commission pursuant to the first subparagraph of Article 5 of Regulation (EU) No 528/2012 of the European Parliament and the Council, in accordance with the criteria that are relevant to human health amongst the criteria established therein: lead (Pb), cobalt (Co). •This product can be safely used in the treatment of children, pregnant or breastfeeding women, or other patient groups considered particularly vulnerable to such substances.

- Do not disassemble or alter the handpiece and cartridge, etc. except as recommended by NSK in this Operation Manual. • Do not drop the handpiece from a height of more than 1 m. • Stop the handpiece from functioning

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